



## KJEFTEKNUSER BB 200

Versatile unit with a focus on fineness

**The Jaw Crusher BB 200 is a floor model designed for a throughput of up to 300 kg/h, suitable for many typical crushing tasks. It accepts feed sizes up to 90 mm. The gap width can be set stepless from 0-30 mm.**

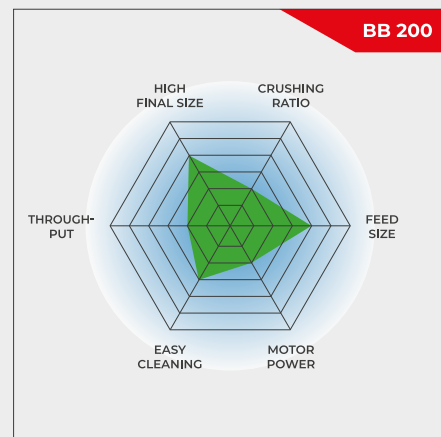
**Depending on the sample material, particle sizes down to 2 mm can be obtained, resulting in a crushing ratio of 45. Its versatile usage includes continuous pre-and fine crushing in one step or process-line versions.**

The RETSCH jaw crusher series is available in 8 different sizes. Throughput and final fineness depend on the crusher type, selected gap width and breaking properties of the sample material.

Robust design, simple handling and cleaning are the features of the RETSCH jaw crusher models. Small amounts of sample are processed batchwise; for larger amounts the crushers can be operated continuously.

## ROBUST & VERSATILE

- | Excellent crushing performance
- | Wide range of materials for contamination free grinding
- | Wear compensation with zero-point adjustment
- | Continuous gap width setting
- | Easy-to-clean crushing chamber
- | Suitable for continuous grinding
- | Connector for dust extraction



Versatile unit with a focus on fineness

KJEFTEKNUSER BB 200

## SUPERIORITY IN DETAIL

Robust design, simple handling and cleaning are the features of the BB 200 jaw crusher. For small sample amounts of up to 5 L the BB 200 can be used batch-wise; for larger amounts it can be operated continuously.

The crushed sample is collected in a removable collector. For larger amounts or continuous crushing operations, the sample collector can be replaced by customer-specific solutions (e.g. a belt conveyor).

A Belleville spring washer integrated in the spindle adjustment provides additional overload protection. The eccentric spindle which moves the crushing arm is driven by a robust brake motor via V-belts. The largest belt pulley also acts as the flywheel to ensure uniform and smooth operation. The BB 200 jaw crusher features maintenance- & lubricant-free sliding bearings.

KJEFTEKNUSER BB 200

## **WEAR COMPENSATION BY ZERO-POINT ADJUSTMENT**

Depending on the material and the throughput, sooner or later the breaking jaws will start to show signs of wear. This means that the set breaking jaw distance or the crushing gap will increase with time. To still be able to obtain reproducible crushing results this wear must be compensated.

RETSCH jaw crushers can be continuously adjusted, allowing for compensation of breaking jaw wear. This is done by slowly altering the gap width setting with the motor running until the breaking jaws are heard to come into contact. The new zero point thus obtained is saved by readjusting the scale.



KJEFTEKNUSER BB 200

## **USER CONVENIENCE COMBINED WITH MAXIMUM SAFETY**

Safety is a top priority with RETSCH jaw crushers. The feed hopper with splash-back protection cannot be accessed by hand. A safety switch and the brake motor ensure an immediate stop if the unit is opened or switched on incorrectly.



For easy cleaning of the crushing chamber, the hinged hopper can be removed in a few simple steps. The jaw crushers run very smoothly and quietly and are virtually maintenance-free.

FOR SAFE AND EFFICIENT CRUSHING PROCESSES

## **ACCESSORIES FOR THE JAW CRUSHER BB 200**

Breaking jaws and wearing plates are available in five different materials to be selected depending on the sample and the analysis to be carried out:



- | manganstål
- | rustfritt stål
- | wolframkarbid
- | stål for tungmetalfri nedmaling

KJEFTEKNUSER BB 200

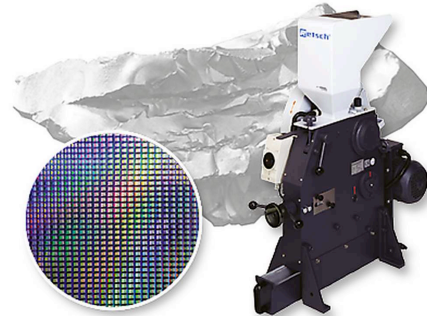
## **JAW CRUSHERS FOR SPECIAL REQUIREMENTS**

Apart from the standard model, the BB 200 jaw crusher is also available in special versions adapted to particular application requirements.

## HEAVY-METAL FREE CRUSHING

All jaw crushers are available in a heavy-metal-free version. The BB 200 can also be supplied in a special version suitable for crushing semiconductor materials. This includes plastic lining of feed hopper and receptacle as well as breaking jaws and wear plates of tungsten carbide.


Thus, the sample does not come into contact with metal materials at any point and no abrasion of the grinding tools impairs the purity of the sample material.



## CONTINUOUS PRE- AND FINE GRINDING

For the rapid, continuous grinding of large quantities of coarse material to analytical fineness, the combination of the jaw crusher BB 200 and the RETSCH Disc Mill DM 200 is the perfect solution.

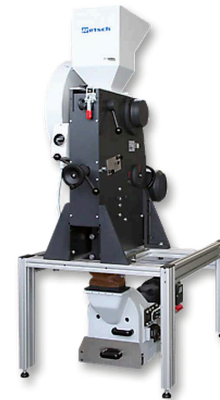The crusher is mounted above the disc mill on a frame and both instruments are connected by a chute. With this construction, samples of up to 90 mm feed size can be ground down to 100 microns in one single step.



## PROCESS-LINE VERSION

The BB 200 jaw crusher is also available in a version suitable for continuous size reduction in online operation, e. g. for quality control during the production process. For these purposes it is supplied without feed hopper and motor protection switch. The voltage of the three-phase AC motor is selected in accordance with customer requirements.



KJEFTEKNUSER BB 200

## TYPICAL SAMPLE MATERIALS

RETSCH's powerful jaw crushers are ideally suited for preliminary crushing of construction materials, ores, granite, oxide ceramics, quartz, slag, silicon, coal, tungsten alloys, cement clinker, etc.



*betong*



*jade*



*glass*



*coal*

To find the best solution for your sample preparation task, visit our application database.

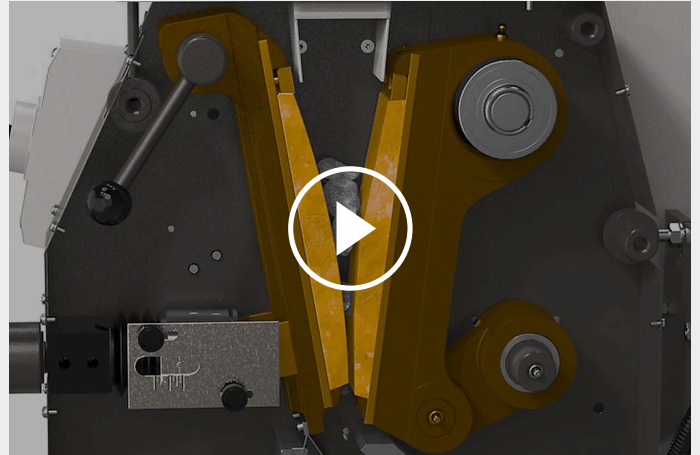
## KJEFTEKNUSER BB 200

### VIRKEMÅTE

Kjefteknuseren BB 200 er en robust og effektiv kjefteknuser med enkelt svinghjul. Prøven passerer gjennom den tilbakeslagssikre trakten og går ned i knusekammeret.

Malingen skjer i den kileformede sjakten mellom den faste og den bevegelige knusearmen, som styres av det eksentriske staget. Den eliptiske bevegelsen gjør at materialet utsettes for trykk og føres nedover.

Når materialet er finere enn den linnstilte spalteåpningen, faller det ned i en uttagbar oppsamlingsbeholder. Den trinnløse spalteinstillingen sikrer en optimal innstilling for å fremstille den ønskede sluttfinheten for prøven.



[Click to view video](#)

KJEFTEKNUSER BB 200

## TECHNICAL DATA

<b>funksjon</b>	grov- og forknusing
<b>bruksområde</b>	byggematerialer, geologi / metallurgi, glass / keramikk, kjemi / plast, maskinbygging / elektroteknikk, miljø / gjenvinning
<b>matemateriale</b>	middels hardt, hardt, sprøtt, seigt
<b>nedmalingsprinsipp</b>	trykk
<b>matestørrelse*</b>	< 90 mm
<b>sluttfinhet*</b>	< 2 mm
<b>Throughput</b>	300 kg/h
<b>material i maleverktøyet</b>	manganstål, rustfritt stål, wolframkarbid, stål 1.1750 (for tungmetallfri nedmaling)
<b>kjeftvidde</b>	90 x 90 mm
<b>spalteåpning</b>	0 - 30 mm
<b>visning av spalteåpning</b>	skala
<b>nullpunktsjustering</b>	ja
<b>vippetrakt</b>	ja
<b>støvsuging:</b>	ja
<b>Maintenance- and lubricant-free sliding bearing</b>	ja
<b>monteringsversjon</b>	valgfritt
<b>oppsamlingsbeholder</b>	5 l
<b>drift</b>	enfasemotor / trefasemotor
<b>drivkraft</b>	1.5 kW
<b>Elektriske tilkoblingsdata</b>	ulike spenninger
<b>Nettilkobling</b>	1-fase / 3-fase
<b>kapslingstype</b>	IP 54
<b>b x h x d lukket</b>	450 x 1160 x 900 mm
<b>nettovekt</b>	~ 300 kg
<b>Normer / standarder</b>	CE

\*avhengig av matemateriale og apparatkonfigurering /apparatinnstillinger

[/bb200](#)

## BESTILLINGSINFO

### JAW CRUSHER BB 200



#### Jaw Crusher BB 200 complete with breaking jaws and wearing plates

				breaking jaws / wearing plates
20.053.0001		BB 200	3/N~ 400 V, 50 Hz	manganese steel / stainless steel
20.053.0002		BB 200	3/N~ 400 V, 50 Hz	stainless steel / stainless steel
20.053.0003		BB 200	3/N~ 400 V, 50 Hz	tungsten carbide / stainless steel
20.053.0043		BB 200	3/N~ 400 V, 50 Hz	tungsten carbide / tungsten carbide
20.059.0001	 	BB 200	3~ 400 V, 50 Hz	tungsten carbide / tungsten carbide, all parts coming into contact with the sample are heavy-metal-free
20.053.0018		BB 200	3/N~ 400 V, 50 Hz	steel 1.1750 / steel 1.0344, for grinding without heavy-metal contamination

other electrical versions available for the same price

### JAW CRUSHER BB 200 VERSION FOR INTEGRATION

Version for integration into a system at the customer's location, with manufacturer's declaration according to EC machine directive, with three-phase motor (voltages on request), without hopper and protective motor switch

20.058.1001		BB 200, manganese steel / stainless steel
20.058.1002		BB 200, stainless steel / stainless steel
20.058.1003		BB 200, tungsten carbide / stainless steel

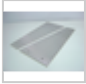

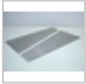
### ACCESSORIES JAW CRUSHER BB 200

1)


#### BREAKING JAWS BB 200

22.048.0004		Breaking jaws manganese steel 1 pair
22.048.0005		Breaking jaws stainless steel 1 pair
22.048.0006		Breaking jaws tungsten carbide 1 pair
22.048.0024		Breaking jaws tungsten carbide 1 pair, for jaw crusher 20.059.0001
22.048.0026		Breaking jaws NiHard 4 1 pair
22.048.0018		Breaking jaws steel 1.1750 1 pair, for grinding without heavy-metal contamination

#### WEARING PLATES BB 200

22.711.0003		Wearing plates stainless steel 1 pair
22.711.0010		Wearing plates tungsten carbide 1 pair
22.711.0006		Wearing plates steel 1.0344 1 pair, for grinding without heavy-metal contamination

#### FURTHER ACCESSORIES BB 200

22.351.0003		V-belt, 1 set (3 pieces)
02.824.0054		Frame for combination BB 200/Disc Mill DM 200